




# CANDELA®



**DERMA LASER**  
Your Laser Service Partner

## Gentle Lase Pro U

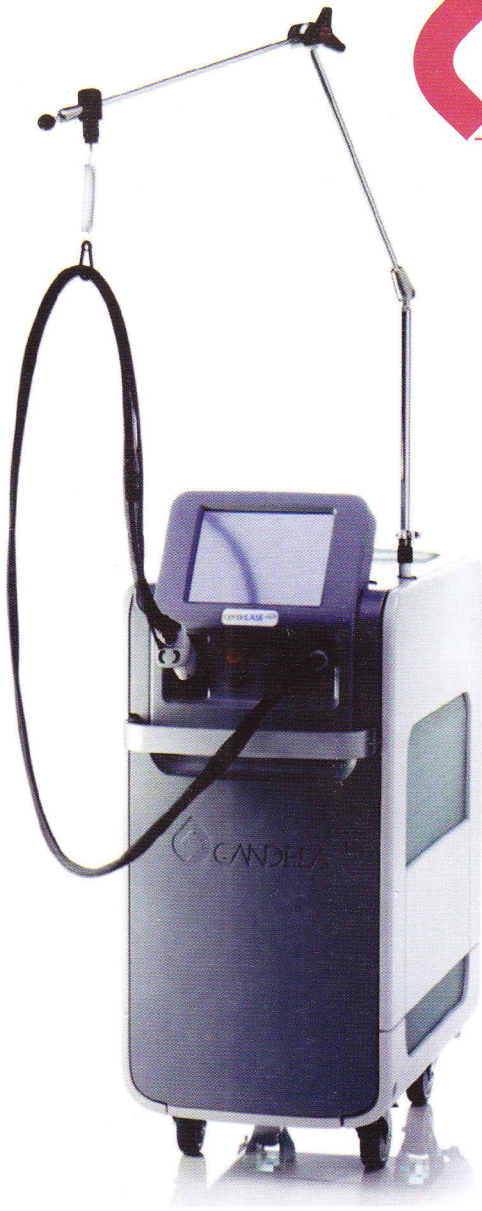

### System Specifications

Laser Type	Alexandrite (Upgrade)
Wavelength	755 nm
Repetition	Up to 10 Hz
Max Delivered Energy	53 joules (J)
Pulse Duration	0.250-100 ms
Spot Sizes	6 mm, 8 mm, 10 mm, 12 mm, 15 mm, 18 mm
Optional Specialty Delivery System Option Spot Sizes	Small - 1.5 mm, 3 mm, 5 mm, 3 x 10 mm Large - 20 mm, 22 mm, 24 mm
Beam Delivery	Lens-coupled optical fiber with handpiece
Pulse Control	Finger switch, foot switch
Dimensions	107 cm H x 46 cm W x 69 cm D (42" x 18" x 27")
Weight	118 kg (260 lbs)
Electrical	200-240 VAC, 50-60 Hz, 20-30 A, 4600 VA single phase (30 A required for highest repetition rates)

#### Faster

An enhanced 2 Hz repetition rate enables you to treat more patients, faster.

#### Choice of Delivery Systems

Choose from multiple ergonomically enhanced hand pieces and lightweight delivery systems to minimize operator fatigue and improve treatment convenience.

#### Intuitive Touch Screen

Select the correct treatment parameters quickly and easily.

#### Variable Pulse Duration

For even greater flexibility in treating vascular and multiple epidermal pigmented lesions.

#### DCD or Air Cooling—It's Your Choice

Choose from either the Air Cooling Compatible™ (ACC™) or the cryogen-based Dynamic Cooling Device™ (DCD™) cooling system, which saves time, is operator independent, scales with fluence and provides consistent protection of the epidermis.

#### Elegant, Smart Design


Sleek, modern lines, a small footprint and lightweight for increased mobility. Deliver fast, effective and comfortable treatments in virtually any environment.

#### Sapphire Protective Windows

Protects system optics, ensuring long lasting performance and the ability to withstand high repetition rates and high energy output.

For Sale & Service Contact 0321-9455593



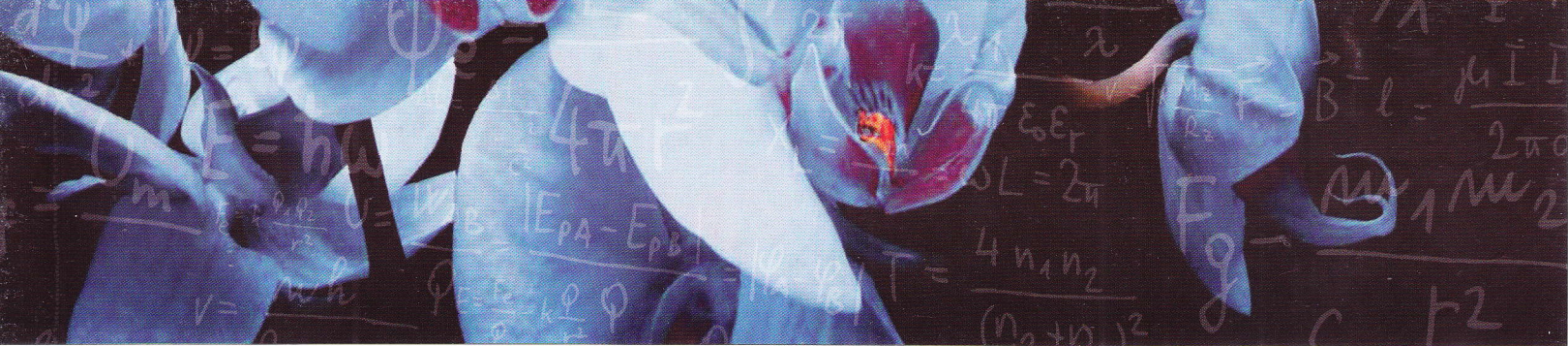

## GentleLase Pro-U Versatility

Practitioners around the world treat more than unwanted hair with GentleLase Pro-U including:

- Hair Removal (ideal for Fitzpatrick skin types I - IV)
- Removal of epidermal pigmented lesions
  - Sun-damaged skin
  - Age spots
  - Freckles
- Treatment of vascular anomalies

## GentleLase Pro-U Upgradable Applications

Expand your GentleLase Pro-U product as your practice expands. Offer the spectrum of hair removal conditions and a wider range of clinical indications:

- Hair Removal on all skin and hair types
- Spider veins
- Telangiectasia and telangiectatic matting on the face and legs

## Clinical Results

### Hair Removal



Before




Post treatment

Photos courtesy of Marcelle C Kutun, M.D., CME, CL, SCLP



Before

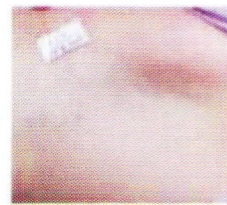

Post treatment

Photos courtesy of Munevver Ercan, M.D.

### Pigmented Lesions



Before



Post treatment

Photos courtesy of Jonathan S. Crance, D.O.

### Vascular Lesions



Before



Post treatment

Photos courtesy of Christine Dierickx, M.D.